

## **TopLine**

Active protection of patients and nursing staff



# **TopLine** – reliable hygiene management for the dirty utility room





Compliance with laws and standards, the duty of care and health and safety in the workplace, regular monitoring and recording of hygiene-related activities: these are just some of the essential requirements placed on the modern nursing profession everywhere in the world. The objective is clear - the prevention of infectious diseases in hospitals and care homes.

MEIKO will provide you with active and professional know-how and innovative, reliable cleaning and disinfection technology: hygiene security from a single source – from the planning stage, technology and chemicals via service right through to professional consultancy in a spirit of partnership, including after-sales assistance.

With our comprehensive hygiene management service which is fine-tuned to the requirements in your premises, we are your dependable partner in the implementation of legal and hygiene-related imperatives.

This is how TopLine will help with the prevention of infections and the protection of patients and staff alike.

- Complete system disinfection at every program start
- Care utensils cooled and dried using HEPA filtered air
- Completely steam-tight system. No discharge of steam at all, even not when the door is opened at the end of the program
- Certified according to EN ISO 15883 by an independent and accredited test institute.
- Energy-optimised disinfection process. Only the energy needed to achieve the pre-set A<sub>0</sub> value is consumed
- Display of the AIM and ACTUAL A<sub>0</sub> value
- A<sub>0</sub> value of 60 and more can be set at will
- Communicative 4 x 20 characters plain text display with energy saving mode
- Nozzle system arranged to provide intensive cleaning, with unique controllable telescopic rotary jet for precisely targeted internal cleaning of the care utensils

# **TopLine** single appliances – the rapid solution for renovation and modernisation projects

Many of the dirty utility rooms of hospitals and care homes no longer comply with current regulations and hygiene standards as well as the demands of nursing staff. TopLine single appliances were specially designed for rapid and uncomplicated installation to

**TopLine 10** – For the clear-of-the-floor wall mounting (with or without slop sink)



The easy way of installing a cleaning and disinfection appliance. Space-saving and ready for connection.

**Dimensions (mm):** w 500, d 450 or 588, h (appliance) 1630



A practical and compact work unit consisting of the TopLine 10 cleaning and disinfection appliance and a slop sink specially designed to accompany it.

Dimensions (mm):

w 1000, d 450 or 588, h (appliance) 1630



replace out-dated cleaning and disinfection technology. They can be modified to fit into the position of the old appliance whatever the make of the original – and to use the same connections.

#### TopLine 20 – Floor mounted appliances (with or without slop sink)



Cleaning and disinfection appliance with plinth. Just place in position and connect. The ideal way of utilizing existing connections.

**Dimensions (mm):** w 500, d 450 or 588, h (appliance) 1730



A work unit which is versatile in its functions, easy to install and compact. Consisting of the TopLine 20 cleaning and disinfection appliance with plinth and a slop sink specially designed to accompany it.

**Dimensions (mm):** w 1000, d 450 or 588, h (appliance) 1730

# TopLine single appliances – seamless and adaptable hygiene security for every work room

#### **TopLine 30** – For the flush wall installation into bathrooms





Cleaning and disinfection appliance TopLine 30 for the bathroom of the patients' rooms, e. g. on the isolation ward. With mounting frame.

Dimensions (mm): Basic frame: w 650, d 450, h 2210 Visible fronts: w 530, h 1760 Cleaning and disinfection appliance TopLine 30 with WC outlet for the bathroom in the patients' room, e. g. on the isolation ward. With mounting frame and complete piping as well as central connections for all commercially available WC types.

Dimensions (mm): Basic frame: w 650, d 450, h 2210 Visible fronts: w 530, h 1760



#### **TopLine 40** – The practical table top model



The TopLine 40 cleaning and disinfection appliance is the perfect complement to the existing cabinet system.

The appliance has a large, uninterrupted work and storage surface.

Upstand available for an additional charge.

**Dimensions (mm): Appliance:** w 900, d 600, h 900



The TopLine 40 cleaning and disinfection appliance with built-in hand wash basin and small storage surface is a suitable and small versatile combined work unit. Upstand available for an additional charge.

Dimensions (mm): Appliance: w 900, d 600, h 1000 Hand wash basin: w 340, d 370, h 150

#### TopLine 40 E - The fastest way to upgrade

Are you renovating or upgrading your facility? The TopLine 40 E is the perfect choice for quickly replacing existing appliances installed under work surfaces. Configurable to fit in all available recesses.

Dimensions (mm): Appliance: w 700, d 588, h 860



# TopLine combined care unit – Hygiene security for all requirements

Hygiene – made to measure. We will work with you to develop an ergonomic utility room design to meet the most demanding requirements of cost-effectiveness, optimisation of work processes, health and safety and hygiene – not to mention any special factors

SAN 14 - TopLine 40 with slop sink / with slop sink and hand wash basin



Combined care unit with work and storage surface, splash screen, slop sink and upstand.

Dimensions (mm): w 1400, d 600, h 900 Slop sink: h 700 Sink: Ø 364



The SAN 14 B compact combined care unit for care needs in the dirty utility room, complete with work and storage surface, splash screen, slop sink, hand wash basin and upstand.

Dimensions (mm): w 1400, d 600, h 1000 Slop sink: h 700 Sink: Ø 364 Hand wash basin: w 340, d 370, h 150



applicable to your premises. We can offer you a hygiene solution customised to meet your required hygiene standards in every room.

#### SAN 19 B - TopLine 40 with slop sink hand wash basin and cabinet

Combined care unit with a large work surface and a cabinet space of 500 mm. With work and storage surface, splash screen, slop sink and upstand.

# Dimensions (mm): w 1900, d 600, h 900 Slop sink: h 700 Sink: Ø 364 Hand wash basin: w 340, d 370, h 150 Cabinet space:

w 500



## SAN 20 BW - TopLine 40 with slop sink, hand wash basin, wash basin and cabinet

Combined care unit with small dimensions and storage surface. With 600 mm cabinet space. With work and storage surface, splash screen, slop sink and upstand.

Dimensions (mm):

w 2000, d 600, h 1000

Slop sink:

h 700

Sink:

Ø 364

Hand wash basin:

w 340, d 370, h 150

Wash basin:

w 500, d 400, h 300

Cabinet space:

w 600



## TopLine combined care units hygiene security for all requirements

#### SAN 24 W - TopLine 40 with slop sink, wash basin, work surface and cabinet 1000 mm

A combined care unit with a large 1000 mm wide cabinet – a very popular choice.

Dimensions (mm): w 2400, d 600, h 900 Slop sink: h 700 Sink: Ø 364 Wash basin: w 500, d 400, h 300 Cabinet space: w 1000



#### SAN 29 BW - TopLine 40 with slop sink, hand wash basin, wash basin and cabinet 1500 mm

Large combined care unit with an ideal length and basin dimensions. All cleaning tasks are optimally combined into a single unit. The basins are ergonomically arranged. The work surface is generously dimensioned.

Dimensions (mm): w 2900, d 600, h 900 Slop sink: h 700 Sink: Ø 364 Hand wash basin: w 340, d 370, h 150 Wash basin: w 500, h 400, d 300 Cabinet space: w 1500



# TopLine ease of operation – Hands-free protection and hygiene




Ease of use and maximum hygiene security for the nursing staff as well as the greatest possible protection for your patients - all guaranteed by the hands-free versions of TopLine equipped with an infrared sensor to operate the door. Also available with an automatic program start on request. The efficient safety stop also protects staff against being accidentally trapped by the door. A foot switch to operate the door is a further option.

## **TopLine 20** & **TopLine 40** – hands-free opening and closing or using the foot switch



TopLine 20 floor model and TopLine 40 table top model with automatic door-opening (optional) and foot switch (optional).

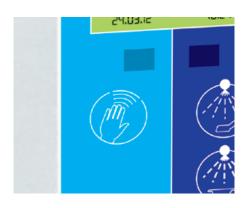

Opening and closing by means of a proximity switch (optional).



Opening and closing by using the foot switch (optional).

### Accessories and range of fixtures



In addition to the MEIKO cleaning and disinfection technology, we can also supply you with a comprehensive range of accessories for safe and ergonomic working and the hygienic storage of care utensils. We illustrate below 6 examples of practical solutions. We would be happy to develop a range of suggested solutions specially for you. Just give us a call. We would be pleased to help you at a moment's notice.



#### Drying and storage racks\*

for disinfected urine bottles, bedpans and receptacles, with storage grid or drip tray. Available in widths of 500 mm to 1.400 mm.



#### Cupboards with doors\*\*

Cupboard units with hinged doors and one intermediate shelf. Available in widths of 500 mm, 600 mm, 800 mm, 1000 mm, 1200 mm and 1500 mm.



#### Wall shelving\*

Rack with height-adjustable shelves available in all shelf lengths. Shelf depth: 350 mm.



#### Cupboards without doors

Cupboards for easy storage with one intermediate shelf. Available in the following widths: 500 mm, 1000 mm and 1200 mm.

- \* Pictures contain special equipment
- \*\* Optional: lamp shade and lamp



#### Cleaning combination\*

Cupboards for installation under work surfaces (complete with drawers and other facilities) can be designed in a range of widths and heights for combined care units.



#### Slop sink\*

Standard model available as an option with integral rim flushing which can be operated independently of the cleaning and disinfection appliance.









720.061.04.08.16/GB/KD/1000

We reserve the right to make changes as part of product development.